


# Mapping for Chronic 64 Channel Probe


ASSY-156 M1

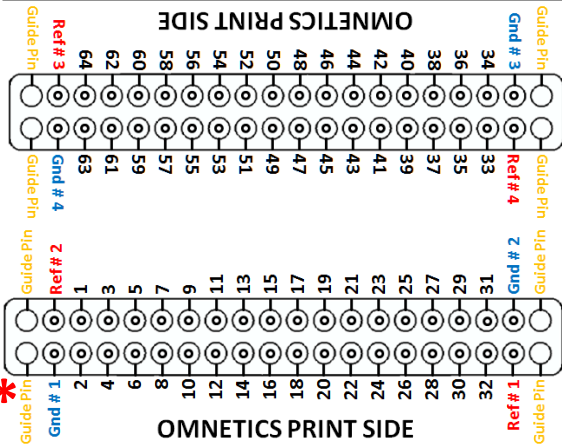

**BACK**

**FRONT**

Red Wire and White Wire correspond to indicated REF and GND pins on the probe connector - verify how these connect to your headstage pre-amplifier (1 wire may be sufficient if you prefer GND / REF shorted-together in your headstage).

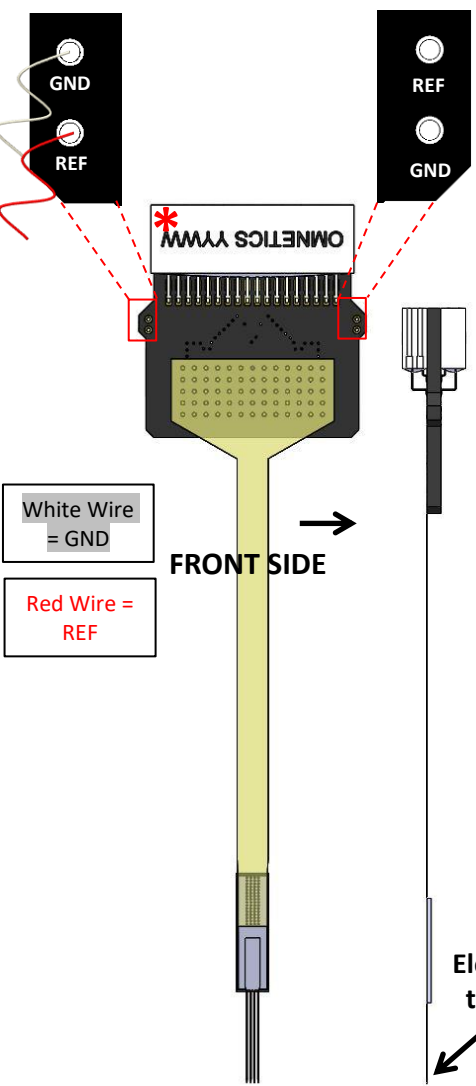

|    |    |    |    |
|----|----|----|----|
| 26 |    | 36 |    |
| 35 | 41 | 38 | 43 |
| 33 | 39 | 25 | 30 |
| 24 | 27 | 34 | 20 |
| 37 | 40 | 21 | 45 |
| 28 | 23 | 22 | 29 |
| 48 | 17 | 47 | 50 |
| 46 | 18 | 19 | 32 |
| 54 | 42 | 56 | 44 |
| 2  | 13 | 16 | 52 |
| 1  | 49 | 15 | 31 |
| 14 | 3  | 4  | 51 |
| 11 | 53 | 64 | 12 |
| 62 | 7  | 9  | 60 |
| 57 | 58 | 55 | 5  |
| 10 | 6  | 63 | 59 |
|    | 61 |    | 8  |

