



# Mapping for Chronic 64 Channel Probe

ASSY-156 P-1 & P-2

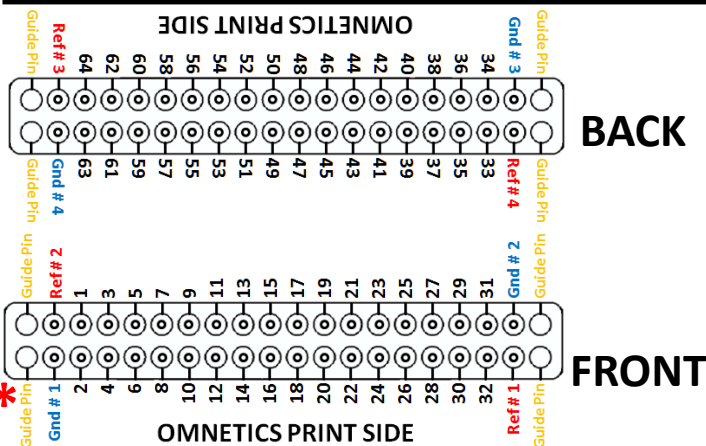

Red Wire and White Wire correspond to indicated REF and GND pins on the probe connector - verify how these connect to your headstage pre-amplifier (1 wire may be sufficient if you prefer GND / REF shorted-together in your headstage).

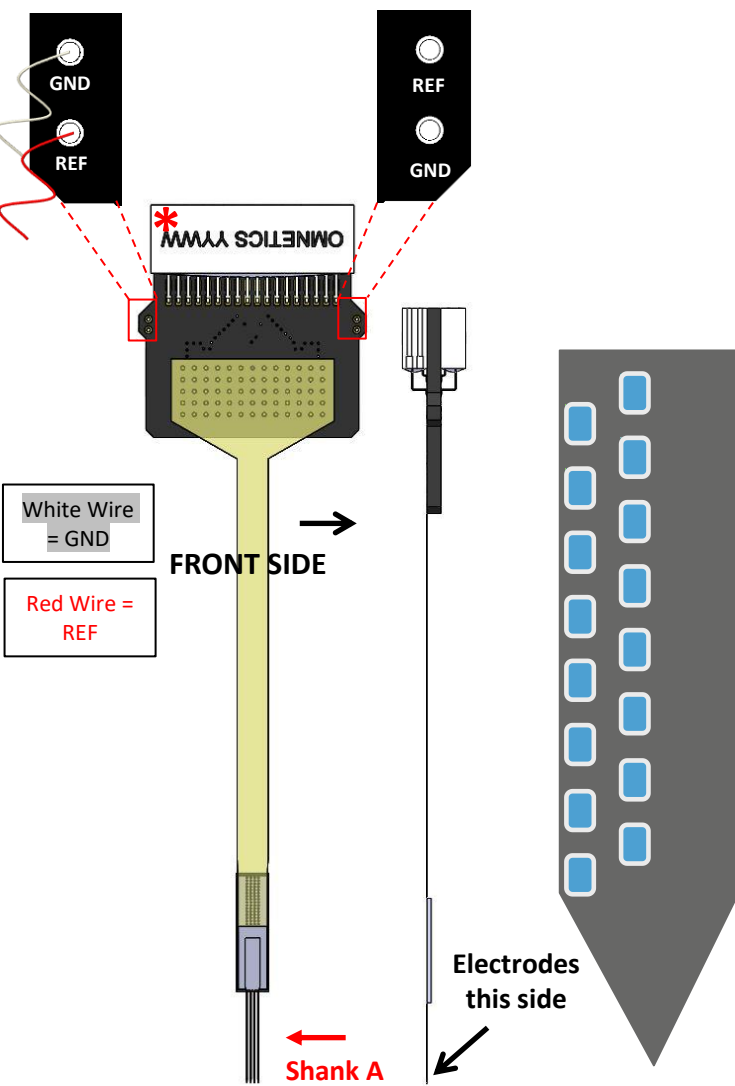

## Shank A

## Shanks B - D

	28		29		3		6
26	41	22	17	12	15	8	56
35	39	43	50	53	1	61	57
33	27	30	48	4	31	63	5
24	40	20	47	14	49	10	58
37	23	45	18	51	16	59	9
36	38	32	46	2	52	62	60
25	34	19	42	13	56	7	64
21		44		54		11	