


Mapping for Y Chronic 64

Channel Probe




ASSY-158 M2 (June 2022 onwards)

Red Wire and White Wire correspond to indicated REF and GND pins on the probe connector - verify how these connect to your headstage pre-amplifier (1 wire may be sufficient if you prefer GND / REF shorted-together in your headstage).

LEFT

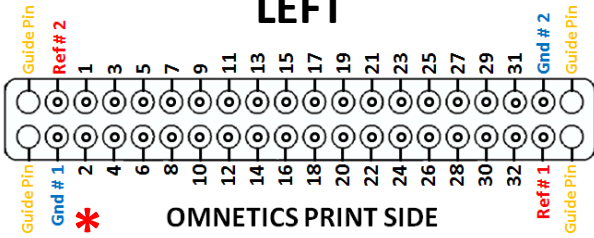

RIGHT

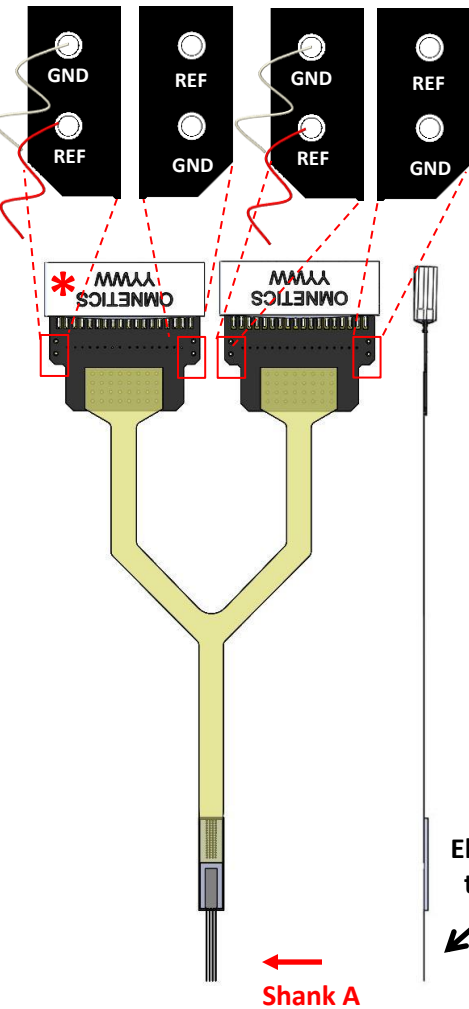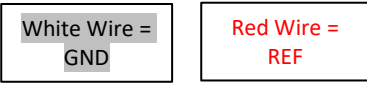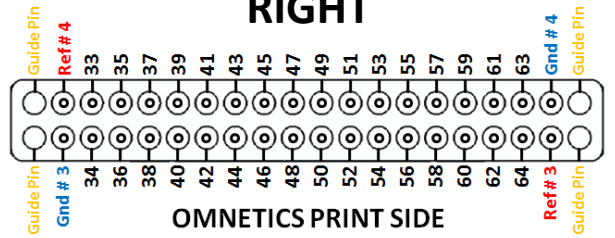

Electrodes
this side

Shank A

| | |
|----|----|
| 64 | 34 |
| 48 | 18 |
| 16 | 54 |
| 61 | 23 |
| 13 | 7 |
| 60 | 38 |
| 29 | 22 |
| 12 | 40 |
| 59 | 56 |
| 32 | 8 |
| 47 | 41 |
| 14 | 55 |
| 31 | 20 |
| 62 | 36 |
| 15 | 6 |
| 63 | 39 |
| 46 | 49 |
| 30 | 19 |
| 58 | 1 |
| 27 | 35 |
| 10 | 4 |
| 42 | 52 |
| 26 | 21 |
| 24 | 5 |
| 25 | 53 |
| 9 | 17 |
| 43 | 33 |
| 57 | 3 |
| 28 | 37 |
| 44 | 51 |
| 11 | 2 |
| 45 | 50 |

Upper-half

Lower-half

