



Mapping for Chronic 64-Channel Probe

Molex connector ASSY-196 H10b

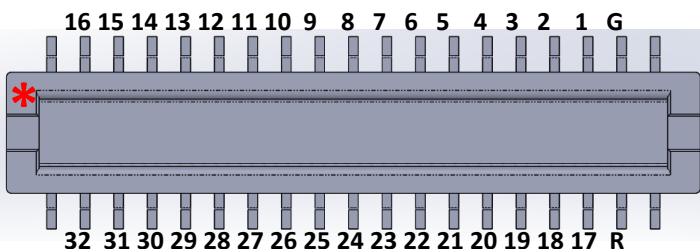

Electrodes are located on the side of the probe as indicated by the arrow below - use Shank A orientation as a reference to construct the channel-map in to your headstage pre-amplifier.

Red Wire and White Wire correspond to indicated REF and GND pins on the probe connector - verify how these connect to your headstage pre-amplifier (often 1 wire is sufficient).

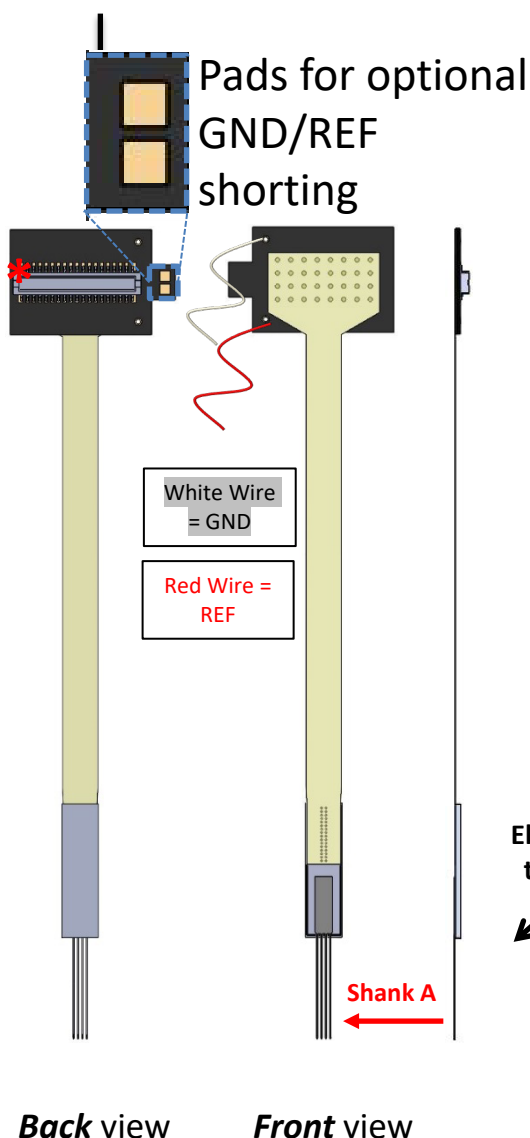

	17	
29		20
	32	
16		1
	2	
31		18
	19	
14		3
	4	
27		22
	21	
12		5
	6	
9		23
	15	
26		8
	30	
10		24
	13	
25		7
	28	
	11	

